# 3696\_TERMOMETRO DIGITAL KELVIN

## FRONT PART

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
## BACK PART

### 3696\_Digital Thermometer Kelvin User manual (Eng)

It is recommended to read this user manual carefully and keep it for future reference.

#### Product Overview
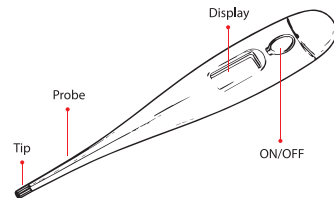
#### Specifications

Measurement range: 32.2 to 43.2 degrees Celsius  
Measuring accuracy: +/2 0.2 degrees Celsius  
Button Battery included LR41

#### Operation

- Press the ON/OFF button. The thermometer will beep once indicating that it is ready to use. The display screen will show the last recorded temperature once and will display LO°C before taking the reading.
- For oral use, place the thermometer under the tongue as indicated below. Keep the mouth closed during reading.

- For an axillary reading, ensure the underarm is dry. Place the thermometer under the arm and bring the arm down against the body.
- The reading will be completed within one minute and the thermometer will beep repeatedly. When the display reading stops flashing and shows consistent temperature, the peak temperature has been reached.
- Read and record the temperature.
- Push the ON/OFF button to turn the unit off. The thermometer will shut off automatically after approximately 10 minutes.

#### Care and storage

- Use the plastic case to store the thermometer.
- Clean the metal tip with isopropyl alcohol or warm soapy water, after each use. Dry thoroughly.
- Wash the probe with warm soapy water after each use, keeping the LCD screen away from the water. Dry thoroughly.
- When the display shows " " or if the digits do not display correctly on the screen, the battery is exhausted and needs replacing. Replace with an LR41 button cell battery.

#### Important information

- Keep out of reach of children.
- Do not allow children to take their own temperature without adult supervision.
- Do not walk or run while taking temperature.
- The body temperature is affected by many factors including exercising, drinking hot or cold drinks prior to taking a reading, etc.
- Oral and axillary temperatures vary by approximately 0.5 degrees Celsius.

#### Warnings

- Do not expose this product to rain or damp conditions.
- Prevent this item from falling as it could cause serious damage.
- Do not dismantle, repair or modify this product by yourself, contact an authorized specialist.
- Do not insert any metal objects into this device as there is a risk of causing a short circuit in the unit.
- Keep out of reach of children.

This symbol both in the product and its user manual indicates that at the end of the devices electrical life, it should be recycled separately from your household waste. There are adequate means to dispose of this material for proper recycling. For more information, contact your local authorities.

Imported by:  
MKTO CATAL IMPORTACIONES S.L. B04267225  
Ctra. De Huerca-Overa, s/n, 04640, Pùbli, Almería, Spain.